

2



DTIC FILE COPY

AD-A203 591

Institute Report No. 315

DTIC
SELECTED
JAN 24 1989
D C

Mutagenic Potential of 1,3-Bis[3-(1-Octoxymethyl)Imidazolium]Propane Dichloride Dihydrate in the Ames *Salmonella*/Mammalian Microsome Mutagenicity Test

Suzanne E. Sebastian, BA, SPC, USA
and
Don W. Korte, Jr., PhD, MAJ, MSC

DISTRIBUTION STATEMENT A

Approved for public release;
Distribution Unlimited

GENETIC TOXICOLOGY BRANCH
DIVISION OF TOXICOLOGY

October 1988

Toxicology Series: 198

LETTERMAN ARMY INSTITUTE OF RESEARCH
PRESIDIO OF SAN FRANCISCO, CALIFORNIA 94129

89 1 23 063

**Mutagenic Potential of 1,3-Bis[3-(1-Octoxymethyl)Imidazolium]Propane Dichloride Dihydrate in the Ames *Salmonella*/Mammalian Microsome Mutagenicity Test (Toxicology Series 198)--
Sebastian and Korte**

This document has been approved for public release and sale; its distribution is unlimited.

Destroy this report when it is no longer needed. Do not return to the originator.

Citation of trade names in this report does not constitute an official endorsement or approval of the use of such items.

This material has been reviewed by Letterman Army Institute of Research and there is no objection to its presentation and/or publication. The opinions or assertions contained herein are the private views of the author(s) and are not to be construed as official or as reflecting the views of the Department of the Army or the Department of Defense. (AR 360-5)



**Edwin S. Beatrice
COL, MC
Commanding**

2 Nov-68
(date)

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE

ADA203591

REPORT DOCUMENTATION PAGE				Form Approved OMB No. 0704-0188	
1a. REPORT SECURITY CLASSIFICATION			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY			3. DISTRIBUTION/AVAILABILITY OF REPORT		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE			Approved for public release; distribution is unlimited.		
4. PERFORMING ORGANIZATION REPORT NUMBER(S) Institute Report No. 315			5. MONITORING ORGANIZATION REPORT NUMBER(S)		
6a. NAME OF PERFORMING ORGANIZATION Genetic Toxicology Branch Division of Toxicology		6b. OFFICE SYMBOL (If applicable) SGRD-ULE-T	7a. NAME OF MONITORING ORGANIZATION Walter Reed Army Institute of Research		
6c. ADDRESS (City, State, and ZIP Code) Letterman Army Institute of Research Presidio of San Francisco, CA 94129-6800			7b. ADDRESS (City, State, and ZIP Code) Washington, D.C. 20307-5100		
8a. NAME OF FUNDING/SPONSORING ORGANIZATION US Army Medical Research & Development Command		8b. OFFICE SYMBOL (If applicable)	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER		
8c. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21701-5012			10. SOURCE OF FUNDING NUMBERS		
			PROGRAM ELEMENT NO. 67234	PROJECT NO. A875	TASK NO. BC
					WORK UNIT ACCESSION NO. DA OH0366
11. TITLE (Include Security Classification) Mutagenic potential of 1,3-bis[3-(1-octoxymethyl)imidazolium]propane dichloride dihydrate in the Ames Salmonella/Mammalian Microsome Mutagenicity Test					
12. PERSONAL AUTHOR(S) Suzanne E. Sebastian, Don W. Korte, Jr.					
13a. TYPE OF REPORT Institute		13b. TIME COVERED FROM 4/21/86 TO 8/23/86		14. DATE OF REPORT (Year, Month, Day) October 1988	
				15. PAGE COUNT 19	
16. SUPPLEMENTARY NOTATION Toxicology Series 198					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)		
FIELD	GROUP	SUB-GROUP	Mutagenicity, Genetic toxicology, Ames Test, 1,3-bis[3-(1-octoxymethyl)imidazolium]propane dichloride dihydrate, oximes, imidazoles, propanes, chlorides. (mamm)		
19. ABSTRACT (Continue on reverse if necessary and identify by block number)					
<p>The mutagenic potential of 1,3-bis[3-(1-octoxymethyl)imidazolium]propane dichloride dihydrate was assessed by using the Ames Salmonella/Mammalian Microsome Mutagenicity Test. Tester strains TA97, TA98, TA100, TA102, TA1537, and TA1538 were exposed to doses ranging from 0.4 mg/plate to 0.000128 mg/plate. The test compound was not mutagenic under conditions of this test.</p>					
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT. <input type="checkbox"/> DTIC USERS			21. ABSTRACT SECURITY CLASSIFICATION Unclassified		
22a. NAME OF RESPONSIBLE INDIVIDUAL Edwin S. Beatrice, COL MC			22b. TELEPHONE (Include Area Code) 415-561-3600		22c. OFFICE SYMBOL SGRD-ULZ

DD Form 1473, JUN 86

Previous editions are obsolete.

SECURITY CLASSIFICATION OF THIS PAGE

UNCLASSIFIED

ABSTRACT

The mutagenic potential of 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE was assessed by using the Ames *Salmonella*/Mammalian Microsome Mutagenicity Test. Tester strains TA97, TA98, TA100, TA102, TA1537, and TA1538 were exposed to doses ranging from 0.4 mg/plate to 0.0000128 mg/plate. The test compound was not mutagenic under conditions of this test.

Key Words: Mutagenicity, Genetic Toxicology, Ames Test, 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE, oxime.



Accession For	
NTIS CRA&I	<input checked="checked" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution	
Availability Codes	
Dist	Avail and/or Special
A-1	

PREFACE

TYPE REPORT: Ames Test GLP Study Report

TESTING FACILITY:

US Army Medical Research and Development Command
Letterman Army Institute of Research
Presidio of San Francisco, CA 94129-6800

SPONSOR:

US Army Medical Research and Development Command
Walter Reed Army Institute of Research
Washington, D.C. 20307-5100

PROJECT/WORK UNIT/APC: 3M162734A875/308/TLEO

GLP STUDY NUMBER: 86002

STUDY DIRECTOR: MAJ Don W. Korte Jr., PhD, MSC

PRINCIPAL INVESTIGATOR: Suzanne E. Sebastian, BA, SPC, USA

REPORT AND DATA MANAGEMENT:

A copy of the final report, study protocol, retired SOP's, stability and purity data on the test compound, and an aliquot of the test compound will be retained in the LAIR Archives.

TEST SUBSTANCE: 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]
PROPANE DICHLORIDE DIHYDRATE

INCLUSIVE STUDY DATES: 21 April 1986 - 23 August 1986

OBJECTIVE:

The objective of this study was to determine the mutagenic potential of 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (LAIR Code TP67) by using the Ames *Salmonella*/Mammalian Microsome Mutagenicity Test.

ACKNOWLEDGMENTS

MAJ John W. Harbell, PhD, MSC; SGT Lillie D. Witcher,
BS; and Ms. Joanne Wong provided research assistance.

SIGNATURES OF PRINCIPAL SCIENTISTS INVOLVED IN THE
STUDY

We, the undersigned, declare that GLP Study 86002 was performed under our supervision, according to the procedures described herein, and that this report is an accurate record of the results obtained.

Don W. Korte, Jr. 27 OCT 88

DON W. KORTE, Jr, PhD / DATE
MAJ, MSC
Study Director

Suzanne E. Sebastian 24 October 88

SUZANNE E. SEBASTIAN, BA / DATE
SPC, USA
Principal Investigator

Conrad Wheeler 27 Oct 88

CONRAD R. WHEELER, PhD / DATE
DAC
Analytical Chemist



DEPARTMENT OF THE ARMY

LETTERMAN ARMY INSTITUTE OF RESEARCH
PRESIDIO OF SAN FRANCISCO, CALIFORNIA 94129-6800

REPLY TO
ATTENTION OF

SGRD-ULZ-QA

1 November 1988

MEMORANDUM FOR RECORD

SUBJECT: GLP Compliance for GLP Study 86002

1. This is to certify that in relation to LAIR GLP Study 86002, the following inspections were made:

15 April 1986	- Protocol Review
21 May 1986	- Plate Incorporation (TP62)
17 March 1987	- Plate Incorporation (TP64)
20 March 1987	- Plate Counting (TP64)

2. The institute report entitled "Mutagenic Potential of 1,3-Bis[3-(1-Octoxymethyl) Imidazolium] Propane Dichloride Dihydrate in the Ames Salmonella/Mammalian Microsome Mutagenicity Test," Toxicology Series 198, was audited on 23 April 1987.

Carolyn M. Lewis

CAROLYN M. LEWIS
Chief, Quality Assurance

TABLE OF CONTENTS

Abstract	i
Preface	iii
Acknowledgments	iv
Signatures of Principal Scientists	v
Report of the Quality Assurance Unit	vi
Table of Contents	vii
INTRODUCTION	1
Objective of the Study	1
MATERIALS AND METHODS	2
Test Compound	2
Test Solvent	2
Chemical Preparation	2
Test Strains	2
Test Format	3
Data Interpretation	4
Deviations from the Protocol/SOP	5
Storage of the Raw Data and Final Report	5
RESULTS	5
DISCUSSION	10
CONCLUSION	10
REFERENCES	11
APPENDICES	12
Appendix A: Chemical Data	13
Appendix B: Individual Plate Scores	14
OFFICIAL DISTRIBUTION LIST	19

**Mutagenic Potential of 1,3-BIS[3-(1-OCTOXYMETHYL)
IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE in the Ames
Salmonella/Mammalian Microsome Mutagenicity Test--
Sebastian and Korte**

INTRODUCTION

1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE was synthesized for a United States Army Medical Research and Development Command program charged with developing more effective oximes for treatment of nerve agent poisoning. The Ames Test is one of a series of tests in which these compounds will be evaluated to determine their relative potential for further development.

The Ames *Salmonella*/Mammalian Microsome Mutagenicity Test is a short-term screening test that utilizes histidine auxotrophic mutant strains of *Salmonella typhimurium* to detect compounds that are potentially mutagenic in mammals. A mammalian microsomal enzyme system is incorporated in the test to increase sensitivity by simulating *in vivo* metabolic activation of the test compound. The Ames Test is an inexpensive yet highly predictive and reliable test for detecting mutagenic activity and thus carcinogenic potential (1).

This evaluation of 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE utilizes a revision of the Ames *Salmonella*/Mammalian Microsome Mutagenicity Test (2). Two new tester strains, a frame-shift strain (TA97) and a strain carrying an ochre mutation on a multicopy plasmid (TA102), are added to the standard tester set.

Objective of the Study

The objective of this study was to determine the mutagenic potential of 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (LAIR Code TP67) by using the revised Ames *Salmonella*/Mammalian Microsome Mutagenicity Test.

MATERIALS AND METHODS

Test Compound

Chemical Name: 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]
PROPANE DICHLORIDE DIHYDRATE

LAIR Code Number: TP67

Physical State: White crystalline solid

Source: SRI International, Menlo Park, CA

Storage: 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE was received from SRI International, 333 Ravenswood Ave., Menlo Park, CA 94025 and assigned the LAIR Code number TP67. The test compound was stored at room temperature (21°C) until used.

Chemical Properties/Analysis: Data provided by SRI International characterizing the chemical composition and purity of the test material are presented in Appendix A along with confirmatory analysis of the test material performed by the Division of Toxicology, LAIR (Presidio of San Francisco, CA).

Test Solvent

The positive control chemicals were dissolved in grade I dimethyl sulfoxide (lot 113F-0450) obtained from Sigma Chemical Co. (St. Louis, MO). The test chemical was dissolved in glass distilled water. Reagent grade water used in this assay was first passed through a Technic Model 301 Reverse Osmosis Unit (Seattle, WA), then through a Corning MP-1 Mega Pure System glass distillation unit (Corning Glass Works, Corning, NY) (3).

Chemical Preparation

On the day of dosing, 100 mg of the test compound was measured into a sterile vial and dissolved in glass distilled water to achieve a 5% (w/v) solution. Aliquots of this solution were used to dose the test plates.

Test Strains

Salmonella strains TA97, TA98, TA100, TA102, TA1537, and TA1538 obtained directly from Dr. Bruce Ames, University of California, Berkeley, were used. These strains were maintained in our laboratory in liquid nitrogen. Quality

control tests were run concurrently with the test substance to establish the validity of their special features and to determine the spontaneous reversion rate. Descriptions of the strains, their genetic markers, and the methods for strain validation are given in the LAIR SOP, OP-STX-1 (4).

Test Format

1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE was evaluated for mutagenic potential according to the revised Ames method (2). A detailed description of the methodology is given in LAIR SOP, OP-STX-1 (4).

Toxicity Tests:

Toxicity tests were conducted to determine a sublethal concentration of the test substance. This toxicity level was found by using minimal glucose agar (MGA) plates, concentrations of 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE ranging from 1.6×10^{-3} mg/plate to 5 mg/plate, and approximately 10^8 cells of TA100 per plate. Top agar containing trace amounts of histidine and biotin was placed on the plates. Strain verification was confirmed on the bacteria, along with a determination of the spontaneous reversion rate. After incubation, the growth on the plates was observed. Since the three highest doses showed a decreased number of macrocolonies (below the spontaneous rate) and an observable reduction in the density of the background lawn, the highest dose selected for the mutagenicity test was 0.04 mg/plate.

Mutagenicity Test:

The test substance was evaluated over a 1000-fold range of concentrations, decreasing from the minimum toxic level (the maximum or limit dose) by a dilution factor of 5, both with and without 0.5 ml of the S-9 microsome fraction. The S-9 (batch R-315) was purchased from Microbiological Associates Inc. (Bethesda, MD). The optimal titer of this S-9, as determined by Microbiological Associates Inc., was 0.75 mg protein/plate. After all the ingredients were added, the top agar was mixed, then overlaid on MGA plates. These plates contained 2% glucose and Vogel Bonner "E" concentrate (5). Plates were incubated upside down in the dark at 37°C for 48 hours. Plates were prepared in triplicate, and the average revertant counts were recorded. The average number of revertants at each dose level was compared to the average number of spontaneous revertants (negative control). The spontaneous reversion rate (with and without S-9) was monitored by averaging the counts from two determinations run

simultaneously with the test compound. The spontaneous reversion rate was determined by inoculating one set of plates before and one set after the test compound plates so that any change in spontaneous reversion rate during the dosing procedure would be detected. This spontaneous reversion rate was also compared with historical values for this laboratory and those cited in Maron and Ames (2). Sterility and strain verification controls were run concurrently. All reagents, test compounds, and media were checked for sterility by plating samples of each on MGA media and incubating them at 37°C with the test plates. The *Salmonella* strains were verified by a standard battery of tests. The integrity of the different *Salmonella* strains used in the assay was verified by the following standard tests:

- Lack of growth (inhibition) in the presence of crystal violet which indicates that the prerequisite alteration of the lipopolysaccharide layer (LP) of the cell wall is present.
- Growth in the presence of ampicillin-impregnated disks which indicates the presence of an ampicillin-resistant R Factor in all strains except TA1537 and TA1538.
- Lack of growth (inhibition) following exposure to ultraviolet light which indicates the absence of the DNA excision-repair mechanism (for all strains except TA102).

Six known mutagens were tested as positive controls to confirm the responsiveness of the strains to the mutation process. Each strain must be tested with at least one positive control but may be tested with several. These compounds, benzo[a]pyrene (lot 79C-05252), 2-aminofluorene (lot 021547), 2-aminoanthracene (lot 020797), mitomycin-C (lot 015F-0655), 4-nitroquinoline-n-oxide (lot 89C-0710) and N-methyl-N'-nitro-N-nitrosoguanidine (lot 127C-0342), were obtained from Sigma Chemical Co. (St. Louis, MO). The test compound and mutagens were handled during this study in accordance with the standards published in NIH Guidelines for the Laboratory Use of Chemical Carcinogens (DHHS Publication No. (NIH) 81-2385, May 1981).

Data Interpretation

According to Brusick (6), a compound is considered mutagenic if a positive dose response (correlated dose response) over three dose concentrations is achieved with at least the highest dose yielding a revertant colony count

greater than or equal to twice the spontaneous colony count for the tester strains TA98 and TA100, or three times the spontaneous colony count for strains TA1537 and TA1538 (2,4). A strong correlated dose response in strain TA100 without a doubling of the individual colony count may also be considered positive.

Maron and Ames (2) consider a compound mutagenic in tester strains TA97 and TA102 if a correlated dose response over three concentrations is achieved with the highest dose yielding a revertant colony count greater than or equal to twice the spontaneous colony count.

Deviations from the Protocol/SOP

Results with TA1535 were not obtained due to technical error. Since the test compound was negative in the other strains, conducting a special assay for TA1535 was deemed unnecessary.

Storage of the Raw Data and Final Report

A copy of the final report, study protocols, raw data, SOPs, and an aliquot of the test compound will be retained in the LAIR Archives.

RESULTS

On 16 May 1986, the toxicity of 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE was determined (Table 1). For this experiment all sterility, strain verification, and negative controls were normal (Table 1). Exposure of the tester strain (TA100) to the three highest doses showed a decrease in the number of macrocolonies, and an observable reduction in the density of the background lawn, indicating chemical toxicity. Therefore, the highest dose selected for the mutagenicity test was 0.04 mg/plate. Normal results were obtained for all sterility and strain verification tests during the Ames Test performed on 20-22 August 1986 (Table 2). Some of the TA100 plates showed excessive clumping of the bacteria, presumably due to inadequate vortexing before plating. An accurate colony count on these plates could not be obtained. The data provided by the remaining strains are sufficient to reach a conclusion on the mutagenic potential of the test compound. 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE did not induce any appreciable increase in the revertant colony counts relative to those of the negative control cultures (Table 3).

TABLE 1: TOXICITY LEVEL DETERMINATION FOR TP67

GLP STUDY NUMBER 86002

<u>TOXICITY DETERMINATION REVERTANT PLATE COUNT (TA100)</u>			
<u>CONCENTRATION</u>	<u>MEAN</u>	<u>±1SD</u>	<u>BACKGROUND LAWN*</u>
START RUN NEGATIVE CONTROL	77	7.5	NL
5.0 mg/plate	0	0	NG
1.0 mg/plate	313	198.8	ST
0.2 mg/plate	206	62.5	ST
0.04 mg/plate	54	4.0	NL
0.008 mg/plate	83	4.6	NL
0.0016 mg/plate	89	10.8	NL
END RUN NEGATIVE CONTROL	92	10.1	NL

STRAIN VERIFICATION FOR TOXICITY DETERMINATION

	<u>TA100*</u>
HISTIDINE REQUIREMENT	NG
AMPICILLIN RESISTANCE	G
UV	NG
CRYSTAL VIOLET SENSITIVITY	NG
STERILITY CONTROL	NG

STERILITY CONTROL FOR TOXICITY DETERMINATION

<u>MATERIAL TESTED</u>	<u>OBSERVATION*</u>
MINIMAL GLUCOSE AGAR PLATES	NG
TOP AGAR	NG
DILUENT WATER	NG
NUTRIENT BROTH	NG
TEST COMPOUND (HIGHEST DOSE)	NG

*NL=Normal Lawn, G=Growth, NG=No Growth, ST=Slight Toxicity

**TABLE 2: STRAIN VERIFICATION AND STERILITY TESTING
FOR THE MUTAGENICITY DETERMINATION OF TP67**

GLP STUDY NUMBER 86002

STRAIN VERIFICATION					
STRAIN	OBSERVATIONS*				
	HISTIDINE REQUIREMENT	AMPICILLIN RESISTANCE	UV REPAIR	CRYSTAL VIOLET	STERILITY CONTROL
TA97	NG	G	NG	NG	NG
TA98	NG	G	NG	NG	NG
TA100	NG	G	NG	NG	NG
TA102	NG	G	G	NG	NG
TA1537	NG	NG	NG	NG	NG
TA1538	NG	NG	NG	NG	NG

STERILITY CONTROL FOR MUTAGENICITY DETERMINATION

<u>MATERIAL TESTED</u>	<u>OBSERVATION*</u>
MINIMAL GLUCOSE AGAR PLATES	NG
TOP AGAR	NG
DILUENT WATER	NG
NUTRIENT BROTH	NG
TEST COMPOUND (HIGHEST DOSE)	NG
S-9	NG

*G = Growth, NG = No Growth

TABLE 3: Mutagenicity Assay for 1,3-BIS[3-(1-OCTOXYMETHYL)
IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (TP67)†

COMPOUND*	DOSE	TA97	TA98	TA100	TA102
WITHOUT S-9					
NEG CONTROL	0.0 mg	96 ±5.8	17 ±3.4	108 ±15.0	76 ±10.4
MITO C	0.5 µg				829 ±117.5
MNNG	2.0 µg			1176 ±50.0	
QNNO	0.5 µg	294 ±14.5			
TP67	40.0 µg	106 ±11.1	23 ±2.5	119 ±23.4	37 ±2.0
TP67	8.0 µg	100 ±8.6	17 ±4.2	91 ±13.5	63 ±12.3
TP67	1.6 µg	98 ±11.5	15 ±0.7		71 ±5.9
TP67	0.32 µg	91 ±3.6	18 ±3.6	114	51 ±7.8
TP67	0.064 µg	114 ±13.1	23 ±4.0	103 ±6.0	71 ±14.0
TP67	0.0128 µg	89 ±6.7	19 ±3.6	109 ±13.6	54 ±11.2
WITH S-9					
NEG CONTROL	0.0 mg	104 ±8.1	26 ±7.5	81 ±9.9	92 ±24.3
AA	2.0 µg		1797 ±51.2	1652 ±28.2	
AF	2.0 µg	264 ±39.7	432 ±12.5		
BP	2.0 µg	333 ±34.6	202 ±27.2		
TP67	40.0 µg	113 ±15.0	24 ±8.3	--	58 ±5.5
TP67	8.0 µg	113 ±16.7	28 ±8.0	90 ±11.3	86 ±3.5
TP67	1.6 µg	110 ±4.6	21 ±0.6	--	93 ±9.9
TP67	0.32 µg	103 ±11.7	21 ±1.0	92 ±9.9	57 ±8.2
TP67	0.064 µg	112 ±11.1	27 ±7.9	105 ±5.7	99 ±10.1
TP67	0.0128 µg	92 ±7.1	28 ±3.8	93 ±15.3	74 ±4.4

† Values represent the mean number of revertants/plate (± standard deviation)

* MITO-C=mitomycin C, MNNG=N-methyl-N'-nitro-N-nitrosoguanidine, QNNO=4-nitroquinoline- n-oxide, AA=2-aminoanthracene, AF=2-aminofluorene, BP=benzo[a]pyrene.

TABLE 3 (cont.): Mutagenicity Assay for 1,3-BIS[3-(1-OCTOXYMETHYL)
IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (TP67)†

COMPOUND*	DOSE/PLATE	TA1537	TA1538
WITHOUT S-9			
NEG CONTROL	0.0 mg	8 ±1.8	14 ±4.6
TP67	40 µg	8 ±2.1	14 ±1.2
TP67	8 µg	6 ±3.1	15 ±5.8
TP67	1.6 µg	7 ±2.6	13 ±1.5
TP67	0.32 µg	6 ±2.1	13 ±4.7
TP67	0.064 µg	7 ±1.2	15 ±1.5
TP67	0.0128 µg	8 ±2.5	12 ±0.6
WITH S-9			
NEG CONTROL	0.0 mg	10 ±3.2	21 ±6.1
AA	2.0 µg	284 ±14.5	1107 ±15.6
AF	2.0 µg		535 ±13.1
BP	2.0 µg	91 ±6.6	105 ±10.7
TP67	40 µg	10 ±0.6	13 ±1.5
TP67	8 µg	11 ±3.8	21 ±3.2
TP67	1.6 µg	8 ±0.6	20 ±6.7
TP67	0.32 µg	14 ±2.1	24 ±3.5
TP67	0.064 µg	19 ±3.8	31 ±6.7
TP67	0.0128 µg	14 ±1.4	26 ±4.0

†Values represent the mean number of revertants/plate (± standard deviation)

*MNNG=N-methyl-N'-nitro-N-nitrosoguanidine, NQNO=4-nitroquinoline-n-oxide, AA=2-aminoanthracene, AF=2-aminofluorene, BP=benzo[a]pyrene.

A tabular presentation of the raw data is included in Appendix B.

DISCUSSION

Certain test criteria must be satisfied before an Ames Test can be considered a valid assessment of a compound's mutagenic potential. First, the special features of the Ames strains must be verified. These features include demonstration of ampicillin resistance, alterations in the LP layer, and deficiency in DNA excision-repair (except TA100). Second, the *Salmonella* strains must be susceptible to mutation by known mutagens. Third, the optimal concentration of the test compound must be determined by treating TA100 with a broad range of doses and observing the potential toxic effects on formation of macrocolonies and microcolonies. If these tests are performed and expected data are obtained, then the results of an Ames Test can be considered valid.

After validation of bacterial strains and selection of optimal sublethal doses, 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE was evaluated in the Ames Test. Criteria for a positive response include both a correlated dose response over three dose concentrations, and a revertant colony count at least two times (TA97, TA98, TA100, TA102) (1,6) or three times (TA1537, TA1538) (2,4) the spontaneous revertant colony count. 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE did not induce the requisite dose-response relationship or the increase in revertant colony counts necessary for a positive response. Thus, the results of this test indicate that 1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE is not mutagenic when evaluated in the Ames test.

CONCLUSION

1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE was evaluated for mutagenic potential in the Ames Test, in both the presence and the absence of metabolic activation, and did not induce a positive mutagenic response under conditions of this study.

REFERENCES

1. Ames BN, McCann J, Yamasaki E. Methods for detection of carcinogens and mutagens with *Salmonella*/Mammalian Microsome Mutagenicity Test. *Mutat Res* 1975;31:347-364.
2. Maron DM, Ames BN. Revised methods for the *Salmonella* Mutagenicity Test. *Mutat Res* 1983;113:173-215.
3. Operation of the Technic Model 301 Reverse Osmosis Pre-Treatment Water System and the Corning Model MP-1 Glass Still. LAIR Standard Operating Procedure OP-STX-94, Presidio of San Francisco, California: Letterman Army Institute of Research, 29 July 1985.
4. Ames *Salmonella*/Mammalian Microsome Mutagenesis Test. LAIR Standard Operating Procedure OP-STX-1, Presidio of San Francisco, California: Letterman Army Institute of Research, 29 August 1986.
5. Vogel HJ, Bonner DM. Acetylornithinase of *E. coli*: Partial purification and some properties. *J Biol Chem* 1956;218:97-106.
6. Brusick D. Genetic toxicology. In: Hayes AW, ed. Principles and methods of toxicology. New York: Raven Press, 1982:223-272.

Sebastian and Korte--12

APPENDICES

APPENDIX A: Chemical Data	13
APPENDIX B: Individual Plate Scores	14


APPENDIX A: Chemical Data

Chemical Name: 1,3-Bis[3-(1-octoxymethyl)imidazolium]propane dichloride dihydrate

SRI Reference Number: 6868-40

LAIR Code Number: TP67

Chemical Structure:

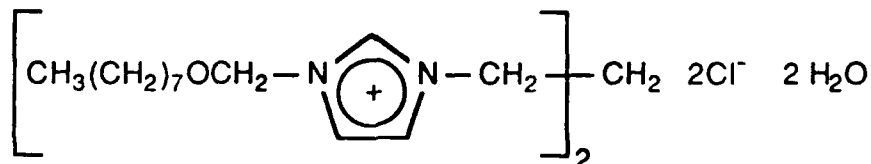

Molecular Formula: $\text{C}_{27}\text{H}_{50}\text{N}_4\text{O}_2\text{Cl}_2 \cdot 2\text{H}_2\text{O}$

Molecular Weight: 569.6

Physical State: White crystalline solid

Analytical Data:

NMR (300 MHz, D_2O): δ 0.92 (t, $J = 6.3$ Hz, 6 H, $\text{CH}_3 - (\text{CH}_2)_7 - \text{O}$), 1.32 (s, 20 H, $\text{CH}_3 - (\text{CH}_2)_5 - \text{CH}_2$), 1.66 (t, $J = 7.2$ Hz, 4 H, $\text{CH}_2 - \text{CH}_2 - \text{O}$), 2.68 (m, $J = 7.5$ Hz, 2 H, $\text{N} - \text{CH}_2 - \text{CH}_2 - \text{CH}_2 - \text{N}$), 3.72 (t, $J = 6.2$ Hz, 4 H, $\text{CH}_2 - \text{CH}_2 - \text{O}$), 4.52 (t, $J = 7.5$ Hz, 4 H, $\text{N} - \text{CH}_2 - \text{CH}_2$), 5.71 (s, 4 H, $\text{O} - \text{CH}_2 - \text{N}$), 7.77 (d, $J = 1.8$ Hz, 2 H, $\text{CH} - \text{CH} - \text{N} - \text{CH}_2 - \text{CH}_2$), 7.80 (d, $J = 1.5$ Hz, 2 H, $\text{OCH}_2\text{NCH} - \text{CH}$).^{*} The NMR spectrum obtained upon receipt of the compound corresponded closely to the spectrum provided by the source (obtained in DMSO). Any discrepancies were due to the difference in solvents as well as the higher field strength and greater resolution of the NMR used to analyze the compound in our lab. No peaks other than those attributable to the compound were observed in the NMR spectrum.

Stability:

NMR data demonstrate that the compound is stable in water (D_2O) for at least 8 days.[†]

Source: Clifford D. Bedford
SRI International
Physical Sciences Division
Menlo Park, CA

^{*}Wheeler CR. Toxicity Testing and Antidotes for Chemical Warfare Agents. Laboratory Notebook #85-12-024, p 13. Letterman Army Institute of Research, Presidio of San Francisco, CA.

[†]*Ibid.* p 1.

APPENDIX B: Individual Plate Scores

1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (TP67)

TOXICITY DETERMINATION WITH TA100

DOSE/PLATE	5.0 mg	1.0 mg	0.2 mg	0.04 mg
PLATE 1	0	325	143	54
PLATE 2	0	505	268	50
PLATE 3	0	108	206	58
background lawn	NG*	ST	NL	NL
DOSE/PLATE	0.008 mg	0.0016 mg	NEG START	NEG END
PLATE 1	84	86	81	90
PLATE 2	78	101	81	103
PLATE 3	87	80	68	83
background lawn	NL	NL	NL	NL

* NG=No Growth, ST=Slight Toxicity, NL=Normal Lawn.

APPENDIX B (cont.): Individual Plate Scores

1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (TP67)

NEGATIVE CONTROL DATA

COMPOUND	DOSE/PLATE	TA97	TA98	TA100	TA102	TA1537	TA1538
<u>WITHOUT S-9</u>							
NEG CONTROL (START RUN)	0.0 mg	93	22	121	61	10	15
		98	17	128	79	8	9
		92	17	89	86	10	18
NEG CONTROL (END RUN)	0.0 mg	98	12	108	66	6	19
		106	14	109	79	9	16
		90	18	94	86	6	8
<u>WITH S-9</u>							
NEG CONTROL (START RUN)	0.0 mg	105	30	*	96	13	19
		111	29	*	108	9	27
		106	20	*	130	9	12
NEG CONTROL (END RUN)	0.0 mg	103	18	*	66	14	29
		88	20	74	72	10	21
		108	37	88	79	5	20

* Plates lost.

APPENDIX B (cont.): Individual Plate Scores

1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (TP67)
POSITIVE CONTROL DATA

COMPOUND	DOSE/PLATE	TA97	TA98	TA100	TA102	TA1537	TA1538
AA	2.0 µg		1812 1839 1740	1631 1684 1641		299 270 284	1098 1125 1098
AF	2.0 µg	219 295 277	422 446 428	* * *	* * *		537 547 521
BP	2.0 µg	369 300 329	216 171 220	* * *		85 98 90	107 114 93
MITO-C	0.5 µg				773 750 964		
MNNG	2.0 µg			1127 1175 1227			
NQNO	0.5 µg	308 279 294					

† AA=2-aminoanthracene, AF=2-aminofluorene, BP=benzo[a]pyrene, MITO-C=mitomycin C,
MNNG=N-methyl-N'-nitro-N-nitrosoguanidine, NQNO=4-nitroquinoline-n-oxide
* Plates lost.

APPENDIX B (cont.): Individual Plate Scores

1,3-BIS(3-(1-OCTOXYMETHYL)IMIDAZOLIUM)PROPANE DICHLORIDE DIHYDRATE (TP67)

MUTAGENICITY DATA WITHOUT S-2

COMPOUND	DOSE/PLATE	TA97	TA98	TA100	TA102	TA1537	TA1538
TP67	40 µg	96	23	146	39	10	13
		118	21	107	37	6	15
		104	26	104	35	7	15
TP67	8 µg	102	14	78	68	9	22
		91	22	105	72	3	12
		108	16	90	49	5	12
TP67	1.6 µg	109	15	*	73	9	13
		98	14	*	75	4	11
		86	*	*	64	8	14
TP67	0.32 µg	95	15	114	56	4	15
		88	17	*	55	7	8
		90	22	*	42	8	17
TP67	0.064 µg	99	25	97	82	8	14
		120	25	103	75	6	15
		123	18	109	55	8	17
TP67	0.0128 µg	82	16	95	66	11	11
		95	18	122	51	8	12
		91	23	111	44	6	12

* Plates lost.

APPENDIX B (cont.): Individual Plate Scores

1,3-BIS[3-(1-OCTOXYMETHYL)IMIDAZOLIUM]PROPANE DICHLORIDE DIHYDRATE (TP67)

MUTAGENICITY DATA WITH S-9

COMPOUND	DOSE/PLATE	TA97	TA98	TA100	TA102	TA1537	TA1538
TP67	40 µg	117	17	*	64	9	12
		125	33	*	55	10	15
		96	21	*	54	10	13
TP67	8 µg	108	36	82	83	8	17
		132	27	98	88	9	23
		100	20	*	*	15	22
TP67	1.6 µg	111	21	*	88	9	18
		105	22	*	104	8	27
		114	21	*	86	8	14
TP67	0.32 µg	94	21	85	64	15	26
		98	20	99	59	12	21
		116	22	*	48	*	*
TP67	0.064 µg	100	30	109	108	22	28
		113	18	101	88	21	27
		122	33	*	101	15	39
TP67	0.0128 µg	98	32	101	79	13	22
		93	25	75	72	*	26
		84	26	100	71	15	30

* Plates lost.

OFFICIAL DISTRIBUTION LIST

Commander
US Army Medical Research
& Development Command
ATTN: SGRD-RMS/Mrs. Madigan
Fort Detrick, MD 21701-5012

Defense Technical Information Center
ATTN: DTIC/DDAB (2 copies)
Cameron Station
Alexandria, VA 22304-6145

Office of Under Secretary of Defense
Research and Engineering
ATTN: R&AT (E&LS), Room 3D129
The Pentagon
Washington, DC 20301-3080

DASG-AAFJML
Army/Air Force Joint Medical Library
Offices of the Surgeons General
5109 Leesburg Pike, Room 670
Falls Church, VA 22041-3258

HQ DA (DASG-ZXA)
WASH DC 20310-2300

Commandant
Academy of Health Sciences
US Army
ATTN: HSHIA-CDM
Fort Sam Houston, TX 78234-6100

Uniformed Services University of
Health Sciences
Office of Grants Management
4301 Jones Bridge Road
Bethesda, MD 20814-4799

US Army Research Office
ATTN: Chemical and Biological
Sciences Division
PO Box 12211
Research Triangle Park, NC 27709-2211

Director
ATTN: SGRD-UWZ-L
Walter Reed Army Institute of Research
Washington, DC. 20307-5100

Commander
US Army Medical Research Institute
of Infectious Diseases
ATTN: SGRD-ULZ-A
Fort Detrick, MD 21701-5011

Commander
US Army Medical Bioengineering Research
and Development Laboratory
ATTN: SGRD-UBG-M
Fort Detrick, Bldg 568
Frederick, MD 21701-5010

Commander
US Army Medical Bioengineering
Research & Development Laboratory
ATTN: Library
Fort Detrick, Bldg 568
Frederick, MD 21701-5010

Commander
US Army Research Institute
of Environmental Medicine
ATTN: SGRD-UE-RSA
Kansas Street
Natick, MA 01760-5007

Commander
US Army Research Institute of
Surgical Research
Fort Sam Houston, TX 78234-6200

Commander
US Army Research Institute of
Chemical Defense
ATTN: SGRD-UV-AJ
Aberdeen Proving Ground, MD 21010-5425

Commander
US Army Aeromedical Research
Laboratory
Fort Rucker, AL 36362-5000

AIR FORCE Office of Scientific
Research (NL)
Building 410, Room A217
Bolling Air Force Base, DC 20332-6448

USAF School of Aerospace Medicine
Document Section
USAFSAM/TSKD
Brooks Air Force Base, TX 78235-5301

Head, Biological Sciences Division
OFFICE OF NAVAL RESEARCH
800 North Quincy Street
Arlington, VA 22217-5000